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OPERATING INSTRUCTIONS

REBOOTS ONE

DE/EN/FR/ES/IT

LIEBER REBOOTS-NUTZER

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ENGLISH

DEAR REBOOTS USER

This guide includes details that should be taken into account by the user, to avoid risks, and to enable the safe use of the device. Should you have any questions about its application, its accessories, or if you have suggestions, please do not hesitate to contact us! We are also happy to help you with any complications that develop during or after using the product.

Telephone: +49 (0) 666 174 891 06

Whatsapp: +49 1573 5990081

Mail: hallo@reboots.de

Web: www.reboots.com

HelpDesk: help.reboots.com

BEFORE YOU GET STARTED

- Please read the operating instructions carefully before using the device.
- Keep this user manual so you can consult it, if necessary.
- Reboots may not be used if you suffer from certain health conditions. Please read the list of contraindications carefully on page 34.

PRODUCT DESCRIPTION


The air pressure massager consists of a device with a pump, cuffs, air pressure sensor, etc., which work together as one unit. The device is connected to the cuffs via a series of hoses. The compression massage is done by cyclically inflating and deflating the air chambers. The pressure in the cuff is controlled by a sensor and microprocessor. The air pressure massager produces sequential compression from distal to proximal, helping to improve the circulation of blood and lymph, prevent DVT and relieve lymphoedema. In cosmetic use, the system provides drainage of loosened fat cells and debris through its massage from the feet towards the upper body.

BEFORE YOUR USE

CONNECT AIR HOSES

Connect the set of hoses to the cuffs, starting with the dark gray connector, and then plugging the other gray connectors onto the appropriate counterparts.

NOTE: We recommend that you do not detach the connectors afterward.

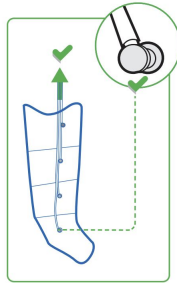DEVICE SET-UP

Ensure that the device is securely placed. Also make sure that the device is free standing and not covered by a blanket, or similar. Do not use the device in damp rooms or those with high humidity, such as in saunas or bathrooms.

TIGHTENING LEG CUFFS

To prevent contamination, we recommend wearing leggings or sweatpants when using the device. Empty your pockets before application!

NOTE: In general, the boots should be pulled up as far as possible. Only the tips of the toes (not the toes!) should stick out at the end. This way you can ensure that your entire legs benefit optimally from the recovery application.

The zippers on the cuffs must be completely closed to avoid the cuffs from suddenly splitting open and causing possible damage.

After plugging it into the power plug, putting on the cuffs, and tightening the zipper, the Reboots application controls on the control unit can be adjusted. You will find more information on the following pages.

NOTE: Always make sure that the hoses are not kinked, and that you are not sitting on the hoses.

ATTACHING THE MULTI-PLUG


Plug the multi-plug into the front of the device as far as it will go. Be sure to insert the multi-plug into the device the right way (the small nipple into the small hole).

REMOTE CONTROL BATTERIES:

Before the initial use: Remove the battery protection and insert the battery (Type: CR2032; Diameter: 20mm, Height: 2.5mm, Voltage: 3V)!

Due to shipping policies, we are not allowed to provide a battery. We kindly ask for your understanding.

WARNING: Keep the remote control and battery out of reach of children as there is a risk of suffocation.

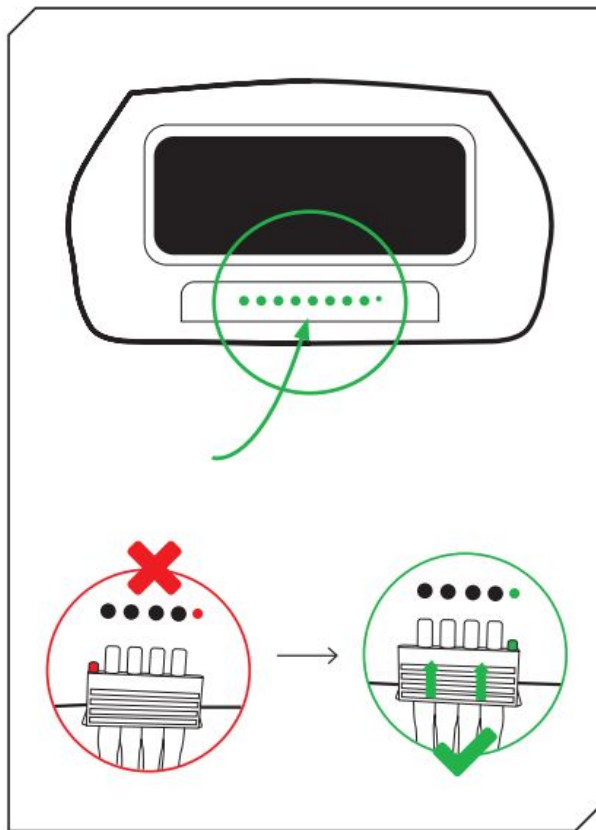
OPERATION

To start a programme, please follow the order of the graphics indicated on the right. During reboot application, hold your legs as horizontally as possible and make sure that the air supply tubes are not pinched. You can adjust the pressure yourself from 30 to 240 mmHg in increments of 30. If necessary, individual chambers can be omitted in order to avoid tender spots, for example.

ATTENTION: Before using the remote control for the first time, the battery protection must be removed and a battery (type: CR2032) must be inserted! Due to shipping regulations, we are unfortunately not allowed to include batteries. We ask for your understanding.

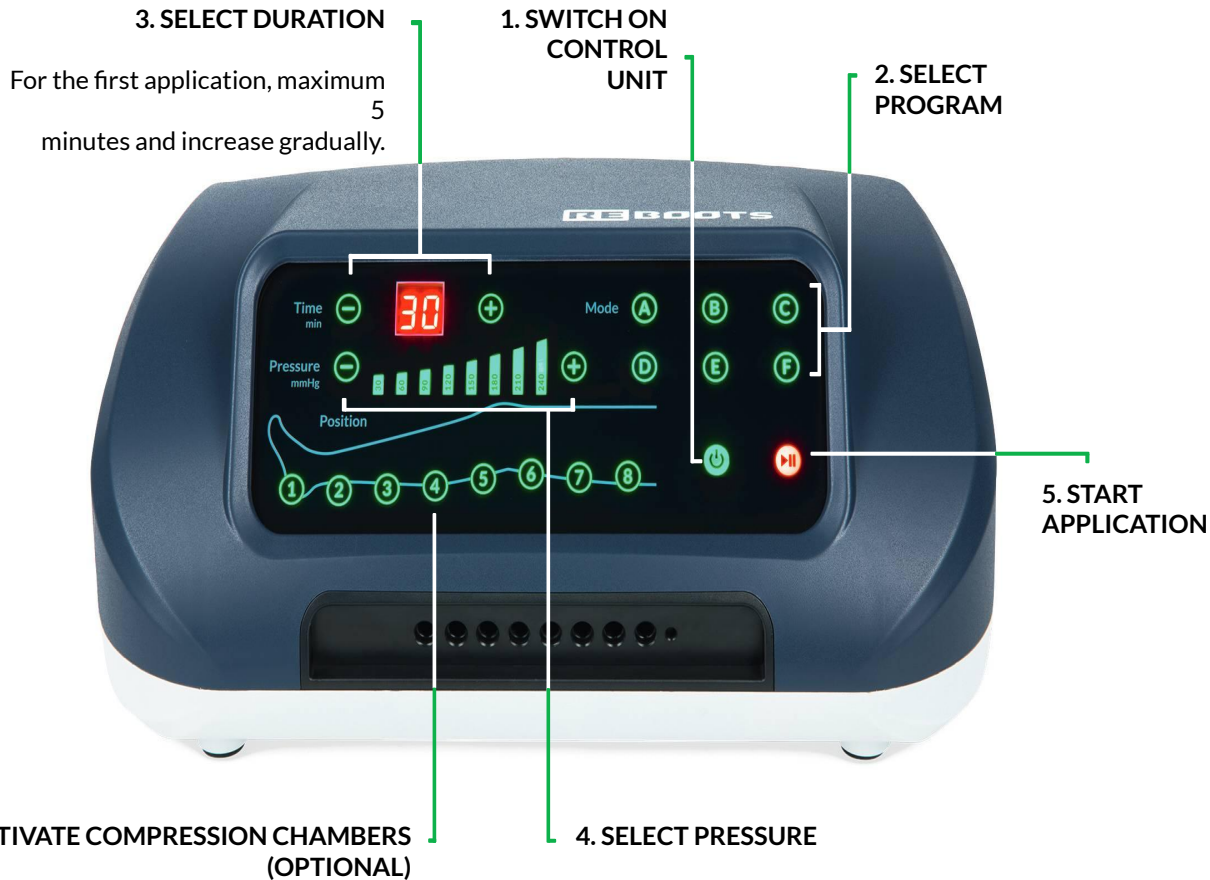Note: Start with low pressure. You can increase during the application!

If you want to exclude certain areas of the leg from the massage, select the compression chambers you want to deactivate here.

PROGRAMME / GUIDE

Maximum regeneration for maximum performance: in order to achieve the best results from your recovery programme, please note the following points:

USE


You can use Reboots both before and after exercising, either to help with relaxation and prevent soreness or to loosen up or muscles in advance. The focus of the massage depends on the programme.

LENGTH AND FREQUENCY OF USE

The optimum length of use after training is between 30 and 50 minutes. Before training, use between 10 and 20 minutes to loosen up your muscles. There is not set maximum length or frequency. In general, you may use your Reboots as long as they feel good. After exercise, it is best to use Reboots on a mid to high level of pressure. Before exercising, you can use Reboots at a mid-level pressure. As long as you feel comfortable during the massage programme, it is fine to use higher pressure or longer programmes. We recommend moving slowly between pressure levels in order to find your optimal setting. More is not always better: the recovery programme should never be unpleasant or painful. Further information about the correct use of this device, programme selection, and other settings can be found online in our Usage Guide. www.reboots.de/blog/reboots-guide


Sequence Example for program

F INTENSE FLOW

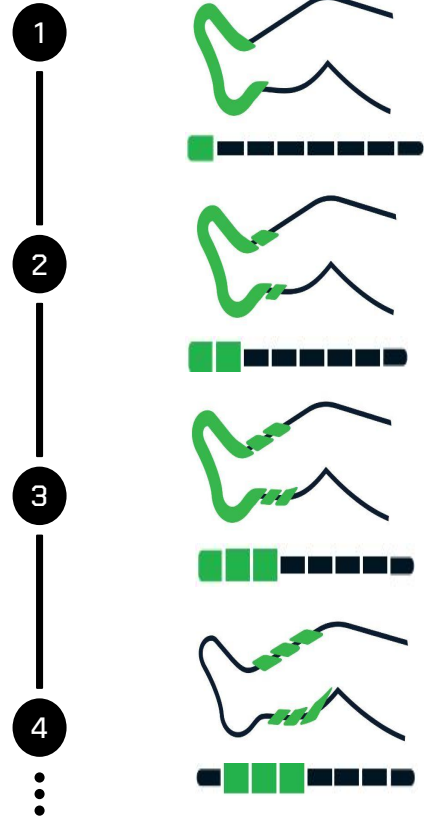
A FLOW

Targeted pressure for muscle relaxation

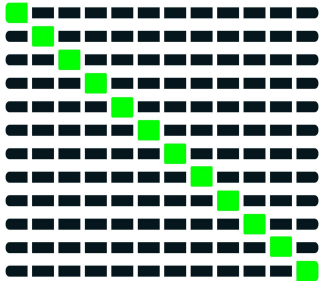
B INTENSE

Continuous pressure for lactate removal

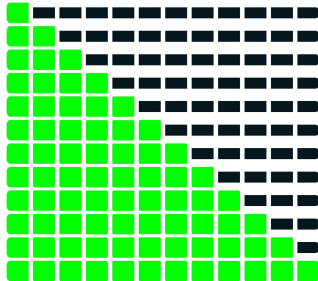
C DOUBLE

Brief applications through more cycles at the same time

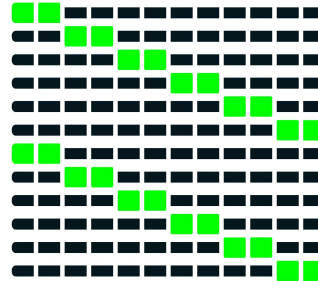
D IMPULSE

Brief applications more intense than program C

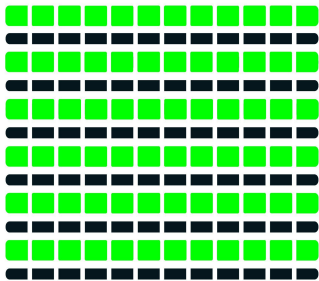
E INTENSE DOUBLE

Wavy massage for lactate removal

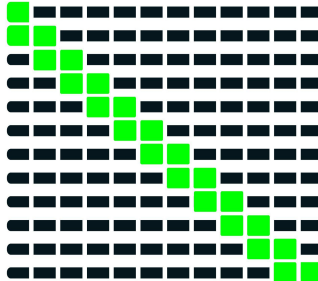F INTENSE FLOW

Con. pressure for lactate removal, less intensive than program B

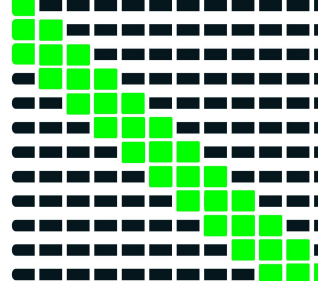AFTER YOUR APPLICATION

COMPLETION & DISMANTLING THE DEVICE

After the set time has elapsed, or after pressing the off switch, you can unplug the multi-plug from the device. If you are not comfortable during the massage, you can also end the massage at any time, as you wish. After that you may open the cuff. To vent the air from the air cushions faster, after application, we recommend pushing the remaining air out of the cuffs with both hands, ensuring that the multi-plug is removed.

After application, we recommend leaving the grey cufflink plugs attached to the cuff and simply removing the multi plug from the device. For storage, place the cuffs on top of each other and loosely roll them together with the attached hose set.

CLEANING

To clean the cuffs and the rest of the accessories, we recommend using a damp cloth. In case of heavy soiling, you can also use a gentle textile cleaner. Before cleaning the cuffs, remove the hose plugs and always allow the cuffs to dry thoroughly. When used by different people, we recommend disinfection of the cuffs, using a textile

disinfectant spray. Use detergents and disinfectants free of oil, benzene, alcohol, gasoline, and/or chemical agents. Do not wash the cuffs in a washing machine or under running water. Switch off the unit before cleaning. Disinfectants such as Descosept pur or comparable products with a similar composition are suitable for disinfecting the cuffs.

STORAGE

Store the product at room temperature (minimum 10°C, maximum 40°C) and in well-ventilated, dry, and dust-free spaces (minimum 30%, maximum 85% relative humidity). Store the product in places free of frost, or it may damage the product.

Also be careful not to damage the cuffs during storage, especially if stored together with sharp objects such as scissors, or objects with sharp edges.

For long-term storage, we recommend that you store the product in the packaging you bought it in.


TROUBLESHOOTING

1. THE DEVICE DOES NOT TURN ON.

If the device does not turn on, please check (A) whether the power plug has been correctly plugged into the socket and/or whether the device is turned on. If the unit is still not turning on, unplug the power plug (B) and unscrew the fuse holder with the inscription “FUSE” on the bottom of the unit (see picture on the left) completely (with a screwdriver or a small coin) and check that the fuse inside is intact. (see picture on the right).

The fuse holder is located on the bottom of the device.

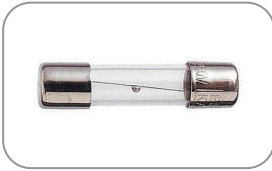Intact fuse: the wire in the fuse is undamaged.


(C) If the fuse is intact, screw the fuse holder with the fuse back into the device completely. If the fuse is damaged, replace the fuse with a commercially-available fine fuse (TR); 3.15 A; 5 x 20 mm). You are also welcome to contact our customer service. (D) Note that the device may only be powered by 220 V power.

2. THE DEVICE IS PUMPING, BUTT ONLY BLOWS ONE OR NONE OF THE TWO CUFFS.

If only one or neither of the two cuffs should inflate after starting the application, please check (A) whether the multi-plug is fully and correctly inserted into the device. Also, make sure (B) that the air hoses are not kinked, for example, because you are sitting on the hoses. Check (C) that the gray cuff plugs are all correctly attached to the cuff.

3. THE CUFFS ARE INFLATING IN THE WRONG ORDER.

Make sure that the multi-plug is correctly inserted into the device. There is a small additional nipple the multi plug, which must be inserted into the hole provided for this purpose. If the multi-plug is inserted into the device the wrong way, the multi-plug will be placed crooked.


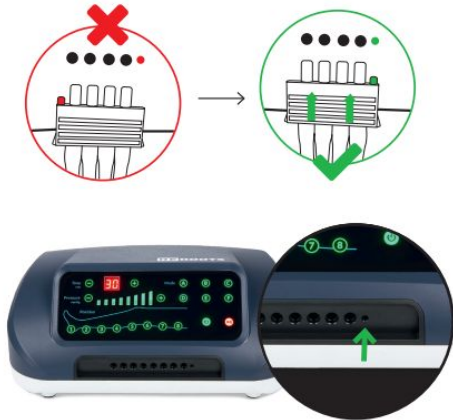
4. MY BOOTS ARE VERY LOUD.

One should be able to talk or relax and chat with another person, without any issues, during the application. If this is not the case, due to the Boots being too noisy, it is an indication that the transport screw on the bottom of the device has not yet been removed. Open the cover on the bottom and remove the safety screw with the enclosed screwdriver. The screw is only used for securing the device during transport, and it can subsequently be disposed of. On short transport routes, we remove the transport screw for you before delivery.



5. YOU CAN HEAR AIR ESCAPING FROM THE DEVICE, HOSES, OR CUFFS.

- (A) Check the hoses and plugs for damage.
- (B) Check that the multi-plug is correctly plugged into the unit.
- (C) Check if the hose is bent or has been pulled off.

INSTRUCTIONS FOR ACCESSORY ITEMS

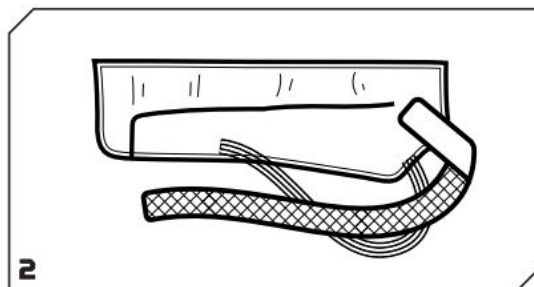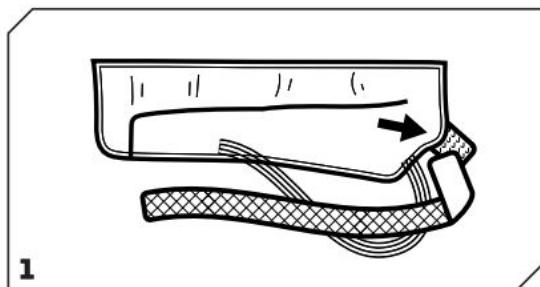Attach the closure flap to the back of the arm cuff.

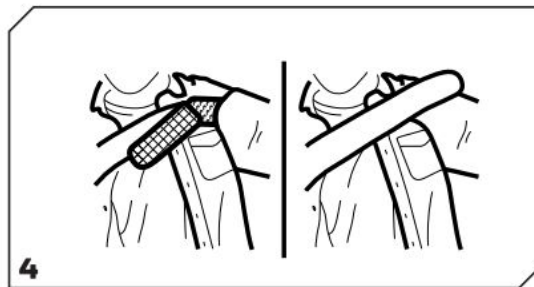
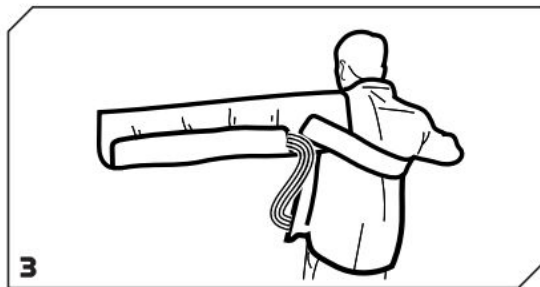
Tighten arm cuff and pay attention to how it fits to the shoulder. Pull the closure flap under the opposite arm/arm pit.

USING INDIVIDUAL CUFFS

Use the dummy plug that sits in the back of the Quick Connector to block unused slots in the Quick Connector.

This allows you to use a single cuff (such as an arm cuff) without having to connect a second cuff to the Quick Connector.

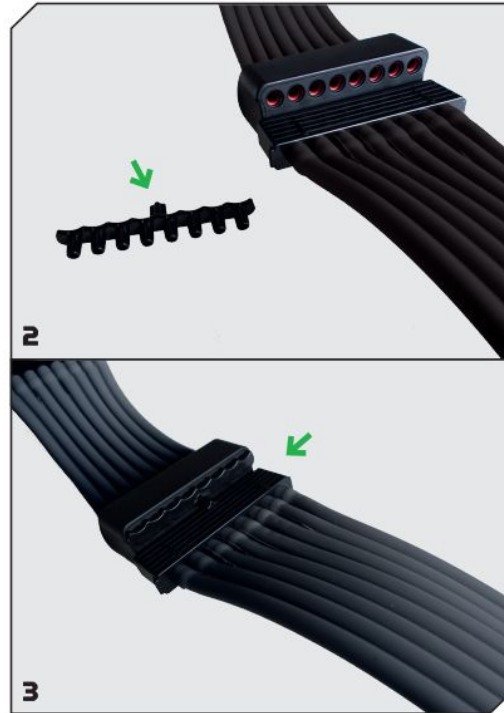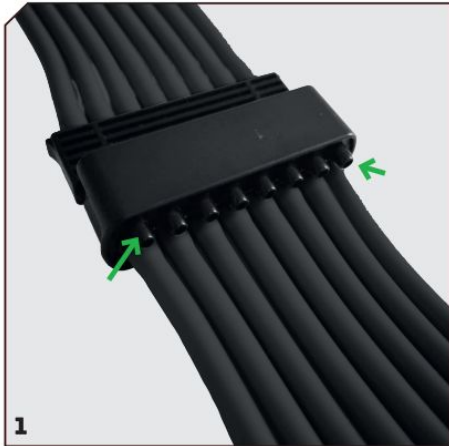INDIVIDUALISING YOUR REBOOTS PANTS 2.0

THE THREE ZIPPERS

Position 1 = outer zipper = narrowest position = for low leg circumferences

Position 2 = middle zipper = middle position = for medium leg circumferences

Position 3 = inner zipper = widest position = for large leg circumferences

There are separate zippers for thighs and lower legs, so you can adjust the trousers very individually.

If none of the possible settings fits exactly, you can create further size gradations with our extensions. Just have a look in our [shop](#).

TECHNICAL DETAILS

Number of chambers One	8
Dimensions control unit [in cm]	30 x 23.7 x 12.6
Weight [in kg]	~ 5
Control way	LED Touch Control mit kabelloser Fernbedienung
Time setting [in min]	5-99
Pressure output [in mmHg]	30-240
Inflation time	10s +/- 1s
Deflation time	50s +/- 2s
Number of programmes	6 Programme
Deactivatable chambers	ja
Mains voltage & frequency [in Hz]	220-230 V, 50-60 Hz

Power output [in VA]	65
Noise level	<65 dB
Protection class	II/ 
Safety class	BF
Waterproof	IP21
Product line	5 Jahre
Application environment	10-40°C, 30-85% (relative), 70kPa-106kPa
Storage environment	-40°C-70°C, 10%-100% (relative), 50kPa-106kPa

REMOTE CONTROL

Battery	DC 3.0V, 1x CR2032
Transmission method	Infrarot
Range	<3m
Size	9,3 x 4,9 x 0,8cm
Weight	20g

PACKAGE CONTENTS

Control unit

Boots / Pants
(M/L/XL)

The scope of delivery of the boots depends on the set purchased

Remote control

Multi-Plug air hose set

DISPOSAL

Note that the product contains small amounts of lead. For this reason, please ensure that the disposal carried out by an official electronic waste disposal service. **IMPORTANT:** Please do not dispose of the product via your general household waste disposal.

ALTERNATIVELY: We are happy to take care of the disposal of your old Reboots free of charge. Just send it to:
Wellcosan GmbH, Hauptstraße 27, 36381 Schlüchtern, Deutschland
Please enclose a note with a request for disposal.

WARRANTY

We offer a 24-month warranty on the control unit. In the event of malfunction, it may be necessary to send the device back for inspection. Please ensure that the shipping box is padded to avoid possible damage during transport. Unfortunately, no claims under the warranty can be brought for defects caused by proper shipping. Terms of warranty: insofar as a legal warranty obligation exists, a replacement will either be provided free of charge or the manufacturer shall have the option to repair the device, under exclusion of the right of cancellation or a reduction in the purchase price. If repeated attempts to effect repair are unsuccessful, or a replacement unit also proves to be defective due to the manufacturer's fault, the customer shall be entitled to cancel the purchase or receive a reduction in the purchase price. Exclusion of warranty: Excluded from the warranty are defects due to damage caused by the use of force, improper operation, externally applied force, or modifications and repair work performed by third parties, e.g. defects caused by incorrectly rated or short-circuited fuses, or defects attributable to normal wear and tear. Conditions of warranty: The warranty claim exists only in connection with the original proof of purchase (invoice). Please keep the proof of purchase in a safe place. Should you have any problems or questions when using the system, please do not hesitate to contact us at +49 (0) 221 958 124 35 or hallo@reboots.de.

LIABILITY DISCLAIMER

Possible contraindications may include: Decompensated heart failure, cardiac output disorders with water in the legs, inflammation of the treated areas of the body, known or suspected leg vein thrombosis, phlebitis and arterial disease, gangrene, dermatitis, arteriosclerosis or venous circulatory disorders, open, untreated or infected wounds and injuries, for leg pain, skin lesions, acute inflammation, malignant tumours, pulmonary oedema, hypertension, extensive thrombophlebitis, thrombosis or suspected thrombosis, erysipelas, severe uncontrolled hypertension, acute soft tissue trauma of the extremities, neuropathy or occlusive processes in the lymphatic drainage system.

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